


CELL LINES

- SARS-CoV-2 (2019-nCoV) Spike D614G Mutant HEK 293 cell line-

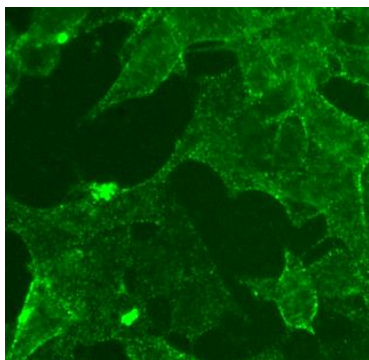


Product Name:	SARS-CoV-2 Spike D614G mutant HEK293 cell line
Catalog Number:	P30910
Cell Line:	HEK293
Resistance:	Hygromycin
Format:	>3x10 ⁶ cells in Cryopreserved vials
Storage:	Liquid Nitrogen

SARS-CoV-2 D614G mutant

spike HEK293

The SARS-CoV-2 spike D614G mutant HEK293 cell line has been developed by stable transfection with SARS-CoV-2 (2019-nCoV) spike D614G mutant protein expression plasmid. SARS-CoV-2 spike HEK293 cell line provides consistent levels of expression of SARS-CoV-2 (2019-nCoV) spike D614G mutant protein in cells surface.

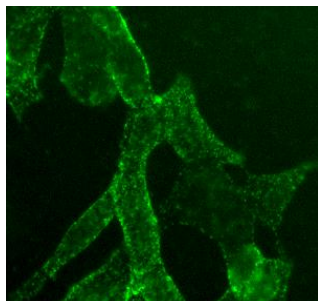

This cell line is intended to be used as an “in vitro” model for research studies.

About SARS-CoV-2 spike protein

The coronavirus S-protein is the structural protein responsible for the crown-like shape of the CoV viral particles.

The SARS-CoV-2 spike protein mediates the membrane fusion process, and utilizes human angiotensin-converting enzyme 2 (hACE2) as the receptor to infect human cells.

Bibliography: Coutard, B., Valle, C., de Lamballerie, X., Canard, B., Seidah, N. G., & Decroly, E. (2020). The spike glycoprotein of the new coronavirus 2019-nCoV contains a furin-like cleavage site absent in CoV of the same clade. *Antiviral research*, 176, 104742.

<https://doi.org/10.1016/j.antiviral.2020.104742>

Xia, S., Liu, M., Wang, C. *et al.* Inhibition of SARS-CoV-2 (previously 2019-nCoV) infection by a highly potent pan-coronavirus fusion inhibitor targeting its spike protein that harbors a high capacity to mediate membrane fusion. *Cell Res* 30, 343–355 (2020). <https://doi.org/10.1038/s41422-020-0305-x>

RT-PCR analysis

The presence of SARS-CoV-2 spike D614G mutant protein mRNA was analyzed by RT-PCR.

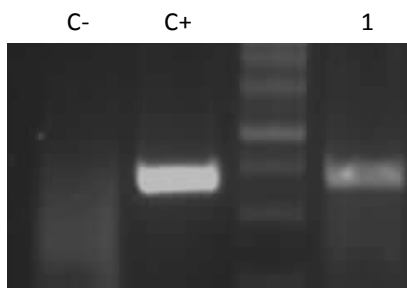


Figure 1. SARS-CoV-2 spike protein RT-PCR analysis. (1) SARS-CoV-2 spike D614G protein HEK293 cell line. Positive Control (C+): SARS-CoV-2 spike protein cDNA. Negative Control (C-): not transfected HEK293 cells.

Immunofluorescence analysis

The detection of SARS-CoV-2 spike D614G protein in the cells surface was carried out by immunofluorescence analysis with a FITC tagged anti-SARS-CoV-2 spike protein antibody.

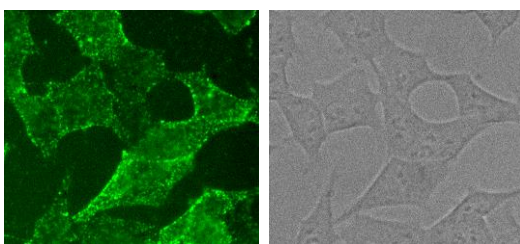


Figure 2. Immunofluorescence assay. The image in the left panel shows the membrane localization of SARS-CoV-2 spike D614G mutant in HEK293 cell line. The image in the right panel shows bright field.

Quality Control

All cells are performance assayed and test negative for mycoplasma, bacteria, yeast and fungi. Cell viability, morphology and proliferative capacity are measured after recovery from cryopreservation. Innoprot guarantees stable expression for many generations and provides support for cell culture and visualization.

THIS PRODUCT IS FOR RESEARCH PURPOSES ONLY. It is not to be used for drug or diagnostic purposes, nor is it intended for human use. Innoprot products may not be resold, modified for resale, or used to manufacture commercial products without written approval of Innovative Technologies in Biological Systems, S.L.